



GoodKnight® H₂O

Effective humidification the whole night
Easy to use
Easy to clean




Improve Comfort and Compliance for CPAP Users

Designed for use with the GoodKnight® family, the GoodKnight® H₂O Heated Humidifier increases air humidity, reducing patient dryness in the nose and throat. By improving comfort, it helps promote patient compliance with positive airway treatment.

In addition, the humidifier features a small footprint that fits neatly with GoodKnight® sleep therapy devices, creating a very small and light CPAP/humidifier combination.

Humid Control Feature

The Humid Control feature ensures a constant temperature difference between the water and ambient air to minimize condensation in the tubing and keep the humidification level constant.

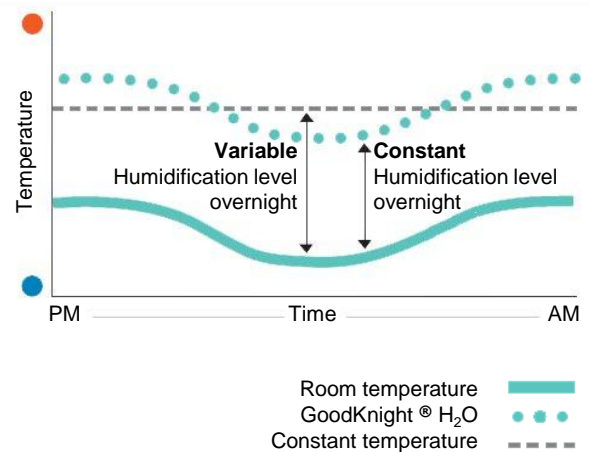

GoodKnight® H₂O is compatible with :
GoodKnight® 420 family
GoodKnight® 425 Bi-Level devices

SPECIFICATIONS :

Performances :

HUMIDIFICATION RATE : 10 mg H₂O/L at 54 L/min and at the maximum setting

TIME OF HEATING : ≤ 45 minutes

SERVICE MAXIMUM PRESSURE : 25 cm H₂O

Physical characteristics :

DIMENSIONS: 170 mm X 95 mm X 280 mm (Including the water chamber)

WEIGHT (WHEN EMPTY) : 0.8 kg

VOLUME USABLE OF THE WATER CHAMBER :
350 ml

MAXIMUM LIQUID LEVEL : Indicated on the chamber

Power characteristics :

INPUT VOLTAGE: Input range 100VAC, 120VAC, 240VAC , 50-60 Hz.

POWER CONSUMPTION : 40 W

IEC 60601-1 CLASSIFICATION : Class II (double insulation)

Environmental requirements :

OPERATING TEMPERATURE : 5°C to 35°C, between 15% and 95% relative humidity (RH), non-condensing

STORAGE TEMPERATURE : -20°C to 60°C, ≤95% relative humidity (RH), non-condensing

ORDERING INFORMATIONS :

Products and accessories :

GoodKnight® H₂O Heated Humidifier

GoodKnight® H₂O water chamber kit

Silicone connector Water Chamber
GoodKnight® H₂O

GoodKnight® H₂O opener water chamber

GoodKnight® H₂O User manual French, German, Italian, Dutch

GoodKnight® 420G with bag and circuit 1.80 m

GoodKnight® 420S with bag and circuit 1.80 m

GoodKnight® 420E with bag and circuit 1.80 m

GoodKnight® 425 with bag and circuit 1.80 m

GoodKnight® 425ST with bag and circuit 1.80 m

References :

M-114600-EU

M-314600-11

M-414600-15

M-414600-12

M-146DFU00-00

M-113900-EU

M-113903-EU

M-113902-EU

M-114500-EU

M-114502-EU



For more information please contact us at the following addresses:
customerservices@sefam-medical.com
marketing@sefam-medical.com

GoodKnight® is a trademark of Nellcor Puritan Bennett Inc.  0459